

**Straumann**

## Iconic tissue level meets immediacy

The Straumann® TLX Implant System is built on a strong basis of scientific studies. It considers key biological principles of hard and soft-tissue healing. It is designed to significantly reduce the risk of inflammation and bone resorption as the implant–abutment interface is moved away from the bone. Straumann® TLX has been perfected for immediacy and is an excellent solution for all other indications to suit the dentist’s preferred treatment protocol—ranging from immediate to conventional placement and loading. The Straumann® TLX Implant System perfectly complements our bone-level BLX Implant System. Both systems use one common drill set and TorcFit™ connection for maximum compatibility with minimum investment.

Simplicity and efficiency

- A one-stage process with restoration at soft-tissue level allows you to use chair time more efficiently,
- Ease of restoration even in the posterior region,
- Highly efficient treatment protocol thanks to straightforward conventional and digitally integrated workflows.

**Institut Straumann AG, Switzerland**

**+41 61 9651111**

**www.straumann.com**

**Dentsply Sirona**

## A new premium solution for an intuitive, reliable and simplified experience

Dentsply Sirona is announcing the launch of the DS OmniTaper Implant System, a new and innovative solution that combines proven technologies with new features. The design is based on the trusted apically tapered Xive Implant design, which has over 20 years of success, modernised with the conical EV connection and OsseoSpeed surface. The DS OmniTaper Implant System will be available in many countries in early 2023.\*

The OmniTaper EV Implant features ActiveBone Control design and has an apically tapered shape. The ActiveBone Control design uses a combination of bone specific preparation protocols and the dual core bone condensing thread design to achieve high primary stability. In combination, this ensures safe and atraumatic implant placement, even in soft bone. The conical EV connection ensures a tight, sta-

ble fit and minimises micromovement and microleakage. In addition, it provides access to the harmonised and comprehensive EV prosthetic range for restorative flexibility. Furthermore, it offers clinicians simplicity and efficiency for immediate chairside solutions.

“We are driven to deliver first-class products and excited to present this new implant system which will further expand clinicians’ ability to meet their patients’ needs,” says Tony Susino, Group Vice President Global Implant Solutions at Dentsply Sirona. “The addition of the OmniTaper EV into the EV family of implants combines the proven implant technologies from Dentsply Sirona into a solution that offers clinicians simplified protocols, efficiency in handling and a portfolio optimised for digital dentistry.”

**Dentsply Sirona, Sweden**

**+46 31 3763000**

**www.dentsplysirona.com/ds-omnitaper**

*\* Due to the different approval and registration times, DS OmniTaper Implant System is not immediately available in all countries.*

EXCELLENCE IN IMMEDIACY

**Straumann® TLX Implant System**  
**Iconic Tissue Level meets**  
**Immediacy.**



Argon Dental

## K3Pro—sustainable tissue retention and impressive prosthetic range

Experienced implantologists have known for a long time: when it comes to the sustainability of hard and soft-tissue stability in two-piece implant systems, the first consideration is the connection between implant and abutment. The so-called Morse taper is particularly solid, a steep and long connection that is equivalent to a cold weld and successfully eliminates any micromovements or gap formations. And thus, reliably ensures tissue preservation. Argon Dental has also relied on this superior concept since the first of its own implants 25 years ago. But with regard to the goal of offering our customers the most modern implant system, the development got a new push in 2010: the self-locking Morse taper became reversible with the help of a threaded channel and an unscrewing instrument that is unique in the industry.



It is not only the prosthetic handling in the practice and laboratory that has benefited enormously—it is above all the options for abutment selection that have resulted from this new development and enable the perfect solution for every indication. In the meantime, you can also perfect your prosthetics completely digitally with K3Pro—for the best that the clever combination of proven principles from the art of engineering and innovative dental prosthetics can make possible for your patients. Celebrate 25 years of K3Pro success with us and be inspired by its unique product features.

**Argon Medical, Germany**  
 +49 6721 3096-0  
[www.argon-medical.com](http://www.argon-medical.com)


BioHorizons Camlog

## Complete solution for bone fixation and membrane stabilisation

Since July 2022, the truFIX system, designed for primary fixation of bone grafting materials, has been distributed by Camlog. The fixation system contains all the necessary components for picking up and screwing in truSCREW (osteosynthesis screw), truTENT (tenting screw) and truTACK (pin) for the fixation of membranes and titanium meshes for reconstruction of bone defects. This comprehensive and user-friendly system makes the use of additional systems and components for defect regeneration unnecessary. TruSCREW with its aggressive self-tapping thread is the ideal bone screw for fixing small bone blocks used in the field of oral and maxillofacial surgery. Thanks to the thread, no pre-drilling is required in most cases. The patented screw design allows effortless insertion into all bone types. The truTENT screw is a complement to truSCREW. The extended shaft and widened head are designed to support a membrane or titanium mesh during augmentation procedures. Smaller, multi-

wall defects can also be successfully regenerated with the help of the tent screws, as they further expand the space and reduce soft-tissue pressure by holding the gingiva over a larger area. Using truTACK, membranes can be stabilised quickly and easily. It has a hexagonal head and threads on the shaft for easy removal. It is inserted like a nail and removed like a screw—a helpful feature.

**Camlog Biotechnologies GmbH, Switzerland**  
 +41 61 5654141  
[www.biohorizonscamlog.com](http://www.biohorizonscamlog.com)

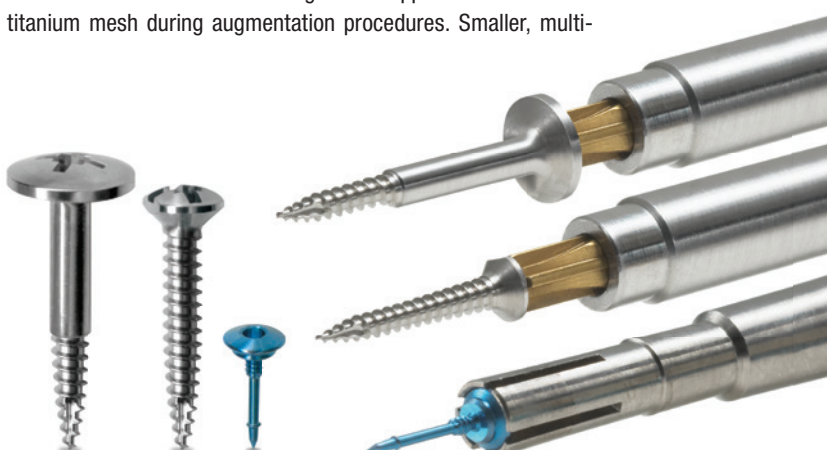

### Titanium meshes

If the augmentation covers more than two tooth widths and has an increasing vertical component, titanium meshes are increasingly being used in addition to the shell technique. They serve as a cage to preserve the space created for the augmentation for regeneration. The titanium meshes are adapted to the defect intra-operatively, filled with augmentation material and fixed in a stable position with screws. They have no barrier function. The titanium meshes are available in different sizes and structures as flat meshes. Depending on the indication, the implantation can be carried out in one or two stages.