


# Use of tilted implants in the treatment of the edentulous posterior maxilla

Author\_Prof Gregor-Georg Zafiropoulos, Germany

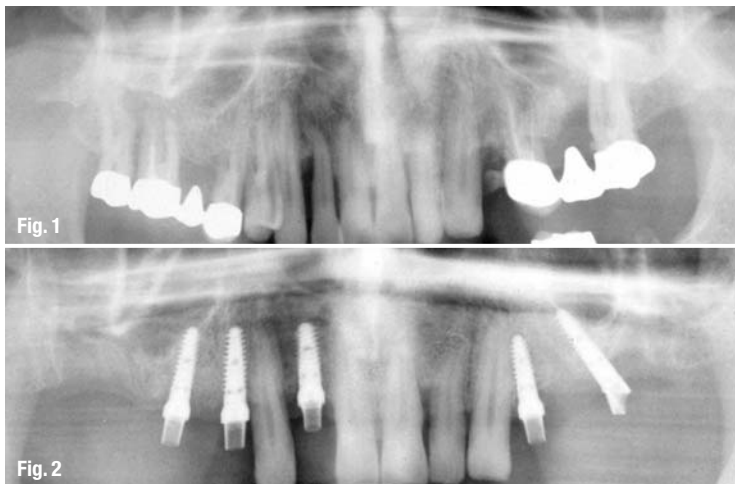

**Fig. 1** \_The orthopantomogram shows the initial situation.

**Fig. 2** \_Individual implant abutments 14 and 16 (both inserted axially).

**Fig. 3** \_Individual implant abutments 24 (inserted axially) and 26 (tilted).

**Fig. 4** \_The orthopantomogram shows implant placement and impression.

\_This case report describes the technique of implant insertion in the edentulous maxillary posterior region at an angle of 35° to avoid a sinus lift procedure and immediate restoration using fixed partial dentures (FPD), e.g. bridges.

Since the introduction of the technique of inserting posterior tilted implants at an angle of up to 35° (in relation to the vertical axis) and the corresponding prefabricated abutments for the treat-

ment of the edentulous maxilla or mandible, implant dentistry has experienced a change in its previously established, conventional surgical and prosthetic thinking (Maló *et al.* 2003, 2005, 2006). With implants inserted in this unorthodox manner, implementation of regenerative measures in the posterior regions of the atrophic and/or partially edentulous maxilla and mandible is avoided, allowing immediate loading and restoration.

Clinical studies show that the success and survival rate of implants inserted at such an angle are comparable to those inserted at a conventional angle (Khatami & Smith 2008; Krekmanov *et al.* 2000; Hinze *et al.* 2010). Further results demonstrated that there was no significant difference in bone loss between implants inserted with a conventional axis and those inserted at an angle (regardless of jaw and/or region; Zampelis *et al.* 2007; Francetti *et al.* 2010).

This implantation technique was developed for the rehabilitation of an edentulous jaw, but only very little information is available about its application for rehabilitation of partial edentulism in posterior regions with FPD (Rocuzzo *et al.* 2009; Cordaro *et al.* 2009). This report presents observa-



tions made over one year of a representative case, in which the edentulous posterior regions of the maxilla were reconstructed using implant-supported FPDs.

## Case report

A 51-year-old male patient (non-smoker) presented himself at the practice of a colleague because of advanced periodontal destruction in the maxillary arch one year before the start of the treatment described in this report (Fig. 1). Teeth #12, 14 to 18, 24 (retained root) and 25 to 27 were extracted. The extraction sockets were covered with dPTFE membranes (Cytoplast, Osteogenics Biomedical) with no additional use of grafting material, as previously described (Hoffmann *et al.* 2008; Zafiroopoulos *et al.* 2010). The mucoperiosteal flap was repositioned and fixed in the region of the papillae using interrupted sutures (Cytoplast, Osteogenics Biomedical). The membranes remained partially exposed and were removed after four weeks. The edentulous areas were then fitted with a model cast prosthesis.

About one year after the extractions in the maxilla, the patient presented at our practice for implant treatment. The patient suffered from bilateral chronic sinusitis and would not allow a sinus augmentation to be performed. Five implants (3.75 mm in diameter and 11.5 mm in length; SoftBone, Dentegris) were placed in regions 12, 14, 16, 24 and 26. The implants inserted in regions 12, 14, 16 and 24 were inserted conventionally, i.e. axially, and an internal sinus lift was performed in region 16 (Fig. 2). The implant inserted in region 26 was inserted at an angle of 35° to the vertical axis and immediately provided with a 35° titanium abutment (DAAS abutment 35°, Dentegris; Fig. 3). An impression was taken using system-specific impression posts (pickup posts were used for the axially placed implants and DAAS posts for the tilted implants, both Dentegris) and a polyether impression material (Impregum, 3M ESPE; Fig. 4). The implants were then provided with system-specific healing caps (Dentegris; Fig. 5).

Three days after implantation, transfer keys were used to fit individual abutments. To fabricate the individual abutments, platinum/iridium/plastic abutments (PTIR abutments, Dentegris) were used as a modelling aid consisting of a prefabricated cast-on base made of platinum-iridium and a screw channel made of residue-free burn-out plastic. To fabricate the abutment for 26, a system-specific castable plastic cylinder (DAAS plastic cylinder, Dentegris) was used. On the same day, both a metal framework made of a cobalt-chromium alloy (ZENOTEC NP, Wieland) and a temporary restoration made of plastic (ZENO-PMMA; Wieland) were milled for immediate restoration of implants 14 to 16 and 24 to 26. The framework was fitted and the temporary restoration was fixed using a temporary cement (TempBond, Kerr; Figs. 6–11). Four months after implant placement and progressive immediate loading with the long-term temporary restoration, the final restoration was fitted using a temporary eugenol-free cement (Implant-Provisional, Alvelogro; Fig. 12).

## Conclusion

Under certain conditions (no active periodontal disease, good patient cooperation, good bone quality), successful early or immediate loading of implants in the posterior maxilla is possible in selected cases. Primary stability and implant design play a major role in the success of the implantation and restoration/loading of bridge restorations on tilted implants (Javed & Romanos 2010, Javed *et al.* 2011). On the basis of the scientific results published to date, it is not possible to discuss evidence-based contraindications (based on the required and measurable values of primary stability, bone density and quality, and influence of occlusal forces). As a result, many questions remain unanswered and risks uncertain (Rocuzzo *et al.* 2009).

In our opinion, the primary benefit of using tilted implants is not necessarily the option of immediate implantation and loading, but firstly the avoidance of augmen-

# TRI

dental implants



**TRI® Performance Concept**  
A revolutionary new implant concept  
in implant dentistry

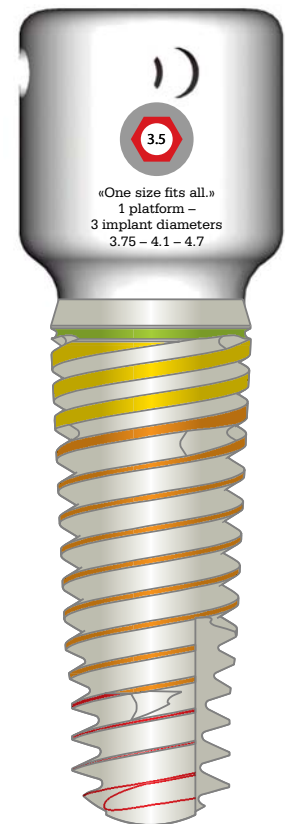


Image:  
Esthetic Healing Collar (TVHC 50-60-A)  
For soft tissue suture fixation, 45 Euro

TRI-Vent Implant (TV41B11), 115 Euro

### 1. TRI® – Friction

This frictional connection guarantees maximum stability and eliminates the risk of a microgap in the platform switch concept.

### 2. TRI® – BoneAdapt

The patented conical implant body has been designed to perform and adapt ideally in the respective bone areas.

### 3. TRI® – Grip

The apical part of this implant allows for a maximum of initial stability for all indications, but specifically for immediate placement.

Enhanced Hard and Soft Tissue Management  
with the patented TRI Dental Implant System.  
Designed in Switzerland. Made in Germany.

### TRI® Performance Days

23.11.2011 Reggio Emilia

26.11.2011 Hamburg

#### Speaker:

Dr. Marius Steigmann

#### Scientific Co-ordinators:

Sergio Spinato

Alberto Caffarri,


Roberto Franceschetti

Test FREE sample!

Through Research Innovative

www.tri-implants.com

Gratis Infoline: 00800 3313 3313

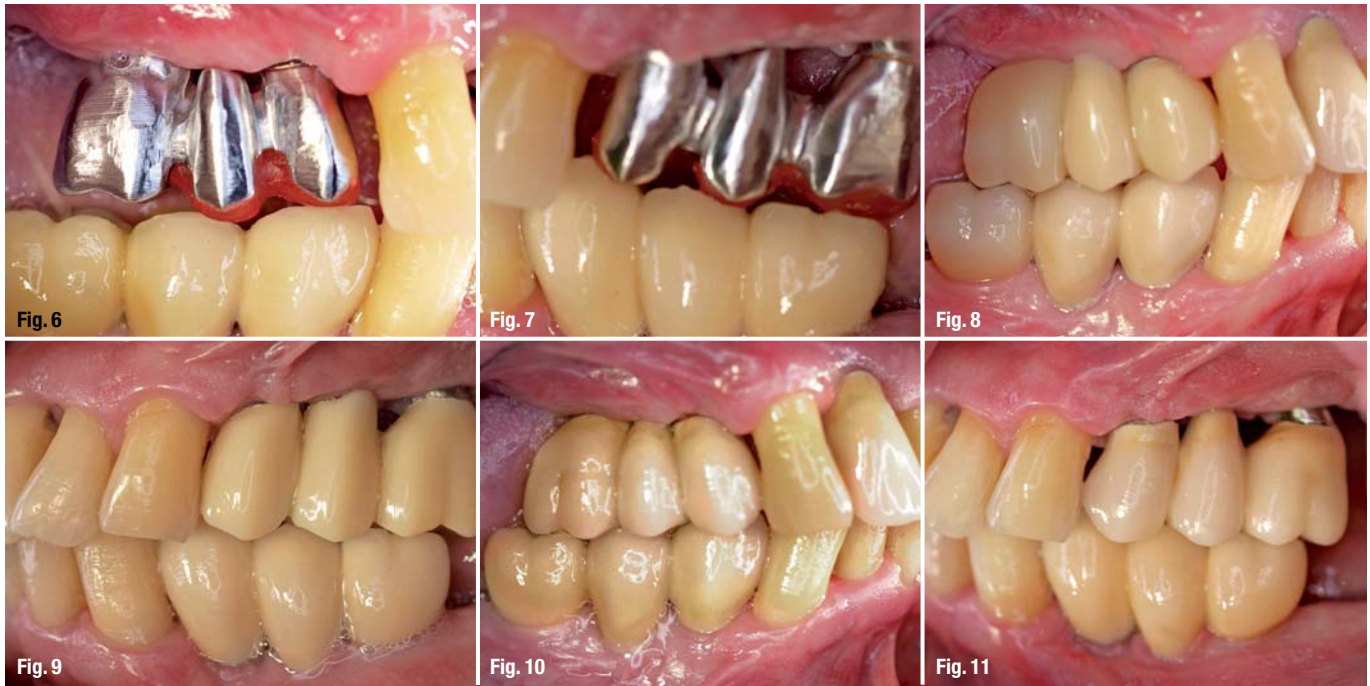

**Fig. 5**\_Implants with healing caps (here an axially placed implant in region 24 and an implant placed at a 35° angle with a DAAS healing cap in region 26).

**Figs. 6 & 7**\_Fitting of the metal frameworks.

**Figs. 8 & 9**\_The implants loaded with temporary FPDs made of ZENO-PMMA.

**Figs. 10 & 11**\_The implants loaded with final metal ceramic FPDs.

tative measures (e.g. sinus lift) and secondly treatment with a fixed restoration (Rosén & Gynther 2007; Aparicio *et al.* 2001). Regardless of the many positive scientific reports, the dentist should be aware of the risk of implant loss associated with immediate implant placement and loading. In addition, the dentist should define his/her own limits and select patients based on strict criteria.

**References**

C. Aparicio, P. Perales & B. Rangert, Tilted implants as an alternative to maxillary sinus grafting: A clinical, radiologic, and periosteal study, *Clinical Implant Dentistry and Related Research*, 3/1 (2001): 39–49.  
 L. Cordaro, F. Torsello & M. Rocuzzo, Implant loading protocols for the partially edentulous posterior mandible, *International Journal of Oral and Maxillofacial Implants*, 24/7 suppl. (2009): 158–168.  
 L. Francetti, D. Romeo, S. Corbella, et al., Bone level changes around axial and tilted implants in full-arch fixed immediate restorations. Interim results of a prospective study, *Clinical Implant Dentistry and Related Research* (published online 26 Oct. 2010) <DOI 10.1111/j.1708-8208.2010.00304.x>

M. Hinze, T. Thalmair, W. Bolz, et al., Immediate loading of fixed provisional prostheses using four implants for the rehabilitation of the edentulous arch: A prospective clinical study, *International Journal of Oral and Maxillofacial Implants*, 25/5 (2010): 1011–8.  
 O. Hoffmann, B.K. Bartee, C. Beaumont, et al., Alveolar bone preservation in extraction sockets using non-resorbable dPTFE membranes: A retrospective non-randomized study, *Journal of Periodontology*, 79/8 (2008): 1355–69.  
 F. Javed & G.E. Romanos, The role of primary stability for successful immediate loading of dental implants. A literature review, *Journal of Dentistry*, 38/8 (2010): 612–20.  
 F. Javed, K. Almas, R. Crespi, et al., Implant surface morphology and primary stability: Is there a connection? *Implant Dentistry*, 20/1 (2011): 40–6.  
 A.H. Khatami & C.R. Smith, "All-on-Four" immediate function concept and clinical report of treatment of an edentulous mandible with a fixed complete denture and milled titanium framework, *Journal of Prosthodontics*, 17/1 (2008): 47–51.  
 L. Krekmanov, M. Kahn, B. Rangert, et al., Tilting of posterior mandibular and maxillary implants for improved prosthesis support, *International Journal of Oral and Maxillofacial Implants*, 15/3 (2000): 405–14.

*Editorial note: The whole list of references is available from the publisher.*



**Fig. 12**\_Orthopantomogram after restoration.

|  |                 |
|--|-----------------|
| <b>_contact</b>  | <b>implants</b> |
| <p><b>Prof Gregor-Georg Zafiropoulos</b><br/>                 Blaues Haus, Sternstr. 61<br/>                 40479 Düsseldorf<br/>                 Germany</p> <p><a href="mailto:zafiropoulos@blaues-haus-duesseldorf.de">zafiropoulos@blaues-haus-duesseldorf.de</a></p> |                 |

”Action is the foundational key  
to all success.“

Pablo Picasso



become an author for “implants”

► Contact Georg Isbaner: [g.isbaner@oemus-media.de](mailto:g.isbaner@oemus-media.de)

