

Endodontic irrigants and irrigant delivery systems

Author_Dr Gary Glassman, Canada

Endodontic treatment is a predictable procedure with high success rates. Success depends on a number of factors, including appropriate instrumentation, successful irrigation and decontamination of the root-canal space to the apices and in areas such as isthmuses. These steps must be followed by complete obturation of the root canals, and placement of a coronal seal, prior to restorative treatment.

Several irrigants and irrigant delivery systems are available, all of which behave differently and have relative advantages and disadvantages. Common root-canal irrigants include sodium hypochlorite (NaOCl), chlorhexidine gluconate, alcohol, hydrogen peroxide and ethylenediaminetetraacetic acid (EDTA). In selecting an irrigant and technique, consideration must be given to their efficacy and safety.

With the introduction of modern techniques, success rates of up to 98% are being achieved.¹ The ultimate goal of endodontic treatment per se is the prevention or treatment of apical periodontitis such that there is complete healing and an absence of infection,² while the overall long-term goal is the placement of a definitive, clinically successful restoration and preservation of the tooth. For these to be achieved, appropriate instrumentation, irrigation, decontamination and root-canal obturation must occur, as well as attainment of a coronal seal. There is evidence that apical periodontitis is a biofilm-induced disease.³ A biofilm is an aggregate of micro-organisms in which cells adhere to each other and/or to a surface. These adherent cells are frequently embedded within a self-produced matrix of extracellular polymeric substance. The presence of micro-organisms embedded in a biofilm and growing in the root-canal system is a key factor for the development of periapical lesions.⁴⁻⁷ Additionally, the root-canal system has a complex anatomy that consists of arborisations, isthmuses and cul-de-sacs that harbour organic tissue and bacterial contaminants (Fig. 1).⁸

The challenge for successful endodontic treatment has always been the removal of vital and necrotic rem-

nants of pulp tissue, debris generated during instrumentation, the dentine smear layer, micro-organisms, and micro-toxins from the root-canal system.⁹

Even with the use of rotary instrumentation, the nickel-titanium instruments currently available only act on the central body of the root canal, resulting in a reliance on irrigation to clean beyond what may be achieved by these instruments.¹⁰ In addition, *Enterococcus faecalis* and *Actinomyces* prevention or treatment of apical periodontitis such as *Actinomyces israelii*—which are both implicated in endodontic infections and in endodontic failure—penetrate deep into dentinal tubules, making their removal through mechanical instrumentation impossible.^{11,12} Finally, *E. faecalis* commonly expresses multidrug resistance,¹³⁻¹⁵ complicating treatment.

Therefore, a suitable irrigant and irrigant delivery system are essential for efficient irrigation and the success of endodontic treatment.¹⁶ Root-canal irrigants must not only be effective for dissolution of the organic of the dental pulp, but also effectively eliminate bacterial contamination and remove the smear layer—the organic and inorganic layer that is created on the wall of the root canal during instrumentation. The ability to deliver irrigants to the root-canal terminus in a safe manner without causing harm to the patient is as important as the efficacy of those irrigants.

Over the years, many irrigating agents have been tried in order to achieve tissue dissolution and bacterial decontamination. The desired attributes of a root-canal irrigant include the ability to dissolve necrotic and pulpal tissue, bacterial decontamination and a broad antimicrobial spectrum, the ability to enter deep into the dentinal tubules, biocompatibility and lack of toxicity, the ability to dissolve inorganic material and remove the smear layer, ease of use, and moderate cost.

As mentioned above, root-canal irrigants currently in use include hydrogen peroxide, NaOCl, EDTA, alcohol and chlorhexidine gluconate. Chlorhexidine

gluconate offers a wide antimicrobial spectrum, the main bacteria associated with endodontic infections (*E. faecalis* and *A. israelii*) are sensitive to it, and it is biocompatible, with no tissue toxicity to the periapical or surrounding tissue.¹⁷ Chlorhexidine gluconate, however, lacks the ability to dissolve necrotic tissue, which limits its usefulness. Hydrogen peroxide as a canal irrigant helps to remove debris by the physical act of irrigation, as well as through effervescing of the solution. However, while an effective anti-bacterial irrigant, hydrogen peroxide does not dissolve necrotic intra-canal tissue and exhibits toxicity to the surrounding tissue. Cases of tissue damage and facial nerve damage have been reported following use of hydrogen peroxide as a root-canal irrigant.¹⁸ Alcohol-based canal irrigants have antimicrobial activity too, but do not dissolve necrotic tissue.

such as NaOCl with EDTA are often used as initial and final rinses to circumvent the shortcomings of a single irrigant.³⁷⁻³⁹ These irrigants must be brought into direct contact with the entire canal-wall surfaces for effective action,^{20,37,40} particularly in the apical portions of small root canals.⁹

The combination of NaOCl and EDTA has been used worldwide for antiseptics of root-canal systems. The concentration of NaOCl used for root-canal irrigation ranges from 2.5 to 6%, depending on the country and local regulations; it has been shown, however, that tissue hydrolysis is greater at the higher end of this range, as demonstrated in a study by Hand et al. comparing 2.5 and 5.25% NaOCl. The higher concentration may also favour superior microbial outcomes.⁴¹ NaOCl has a broad antimicrobial spectrum,²⁰

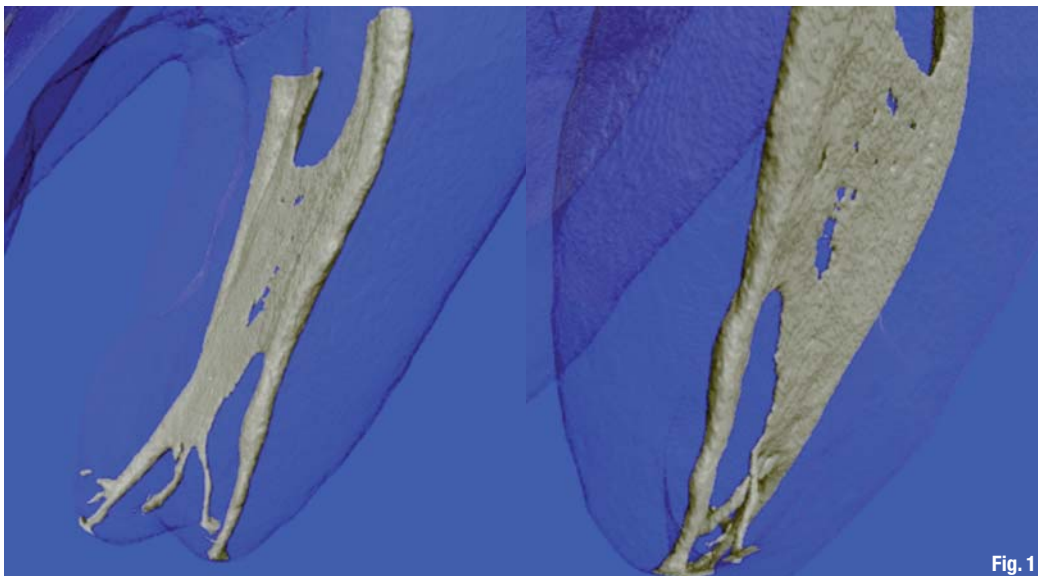


Fig. 1 _ Root-canal complex.

(Image courtesy of Dr Ronald Ordinala Zapata, Brazil.)

<http://www.facebook.com/TheInternalAnatomyOfTheHumanTeeth>

The irrigant that satisfies most of the requirements for a root-canal irrigant is NaOCl.^{19,20} It has the unique ability to dissolve necrotic tissue and the organic components of the smear layer.^{19,21,22} It also kills sessile endodontic pathogens organised in a biofilm.^{23,24} There is no other root-canal irrigant that can meet all these requirements, even with the use of methods such as lowering the pH,²⁵⁻²⁷ increasing the temperature,²⁸⁻³² or adding surfactants to increase the wetting efficacy of the irrigant.^{33,34} However, although NaOCl appears to be the most desirable single endodontic irrigant, it cannot dissolve inorganic dentine particles and thus cannot prevent the formation of a smear layer during instrumentation.³⁵

Calcifications hindering mechanical preparation are frequently encountered in the root-canal system, further complicating treatment. Demineralising agents such as EDTA have therefore been recommended as adjuvants in root-canal therapy.^{20,36} Thus, in contemporary endodontic practice, dual irrigants

including but not limited to *E. faecalis*. NaOCl is superior among irrigating agents that dissolve organic matter. EDTA is a chelating agent that aids in smear layer removal and increases dentine permeability,^{42,43} which will allow further irrigation with NaOCl to penetrate deep into the dentinal tubules.⁴⁴

_General safety precautions

Regardless of which irrigant and irrigation system is employed, and particularly if an irrigant with tissue toxicity is used, there are several general precautions that must be followed. A rubber dam must be used and a good seal obtained to ensure that no irrigant can spill from the pulp chamber into the oral cavity. If deep caries or a fracture is present adjacent to the rubber dam on the tooth being isolated, a temporary sealing material must be used prior to performing the procedure to ensure a good rubber dam seal. It is also important to protect the patient's eyes with safety glasses and protect clothing from irrigant splatter or spill.

It is very important to note that while NaOCl has unique properties that satisfy most requirements for a root-canal irrigant, it also exhibits tissue toxicity that can result in damage to the adjacent tissue, including nerve damage should NaOCl incidents occur during canal irrigation. Furthermore, Salzgeber reported in the 1970s that apical extrusion of an endodontic irrigant routinely occurred *in vivo*.⁴⁵ This highlights the importance of using devices and techniques that minimise or prevent this. NaOCl incidents are discussed later in this article.

_Irrigant delivery systems

Root-canal irrigation systems can be divided into two categories: manual agitation techniques and machine-assisted agitation techniques.⁹ Manual irrigation includes positive-pressure irrigation, which is commonly performed with a syringe and a side-vented needle. Machine-assisted irrigation techniques include sonics and ultrasonics, as well as newer systems such as the EndoVac (SybronEndo), which delivers apical negative-pressure irrigation,⁴⁶ the plastic rotary FFile (Plastic Endo),^{47, 48} the Vibringe (Vibringe),⁴⁹ the Rinsendo (Air Techniques),⁹ and the EndoActivator (DENTSPLY Tulsa Dental Specialties).⁹ Two important factors that should be considered during the process of irrigation are whether the irrigation system can deliver the irrigant to the whole extent of the root-canal system, particularly to the apical third, and whether the irrigant is capable of debriding areas that could not be reached with mechanical instrumentation, such as lateral canals and isthmuses. When evaluating irrigation of the apical third, the phenomenon of apical vapour lock should be considered.⁵⁰⁻⁵²

_Apical vapour lock

Since roots are surrounded by the periodontium, and unless the root-canal foramen is open, the root canal behaves like a close-ended channel. This produces an apical vapour lock that resists displacement during instrumentation and final irrigation, thus preventing the flow of irrigant into the apical region and adequate debridement of the root-canal system.^{53, 54} Apical vapour lock also results in gas entrapment at the apical third.⁹ During irrigation, NaOCl reacts with organic tissue in the root-canal system, and the resulting hydrolysis liberates abundant quantities of ammonia and carbon dioxide.⁵⁵ This gaseous mixture is trapped in the apical region and quickly forms a column of gas into which further fluid penetration is impossible. Extension of instruments into this vapour lock does not reduce or remove the gas bubble,⁵⁶ just as it does not enable adequate flow of irrigant.

The phenomenon of apical vapour lock has been confirmed in studies in which roots were embedded

in a polyvinylsiloxane impression material to restrict fluid flow through the apical foramen, simulating a close-ended channel. The result in these studies was incomplete debridement of the apical part of the canal walls with the use of a positive-pressure syringe delivery technique.⁵⁷⁻⁶⁰ Micro-CT scanning and histological tests conducted by Tay et al. have also confirmed the presence of apical vapour lock.⁶⁰ In fact, studies conducted without ensuring a close-ended channel cannot be regarded as conclusive on the efficacy of irrigants and the irrigant system.⁶¹⁻⁶³ The apical vapour lock may also explain why in a number of studies investigators were unable to demonstrate a clean apical third in sealed root canals.^{59, 64-66}

In a paper published in 1983, based on research Chow determined that traditional positive-pressure irrigation had virtually no effect apical to the orifice of the irrigation needle in a closed root-canal system.⁶⁷ Fluid exchange and debris displacement were minimal. Equally important to his primary findings, Chow set forth an infallible paradigm for endodontic irrigation: "For the solution to be mechanically effective in removing all the particles, it has to: (a) reach the apex; (b) create a current (force); and (c) carry the particles away."⁶⁷ The apical vapour lock and consideration for the patient's safety have always prevented the thorough cleaning of the apical 3 mm. It is critically important to determine which irrigation system will effectively irrigate the apical third, as well as isthmuses and lateral canals,¹⁶ and in a safe manner that prevents the extrusion of irrigant.

_Manual agitation techniques

By far the most common and conventional set of irrigation techniques, manual irrigation involves dispensing of an irrigant into a canal through needles/cannulae of variable gauges, either passively or with agitation by moving the needle up and down the canal space without binding it on the canal walls. This allows good control of needle depth and the volume of irrigant that is flushed through the canal.^{9, 63} However, the closer the needle tip is positioned to the apical tissue, the greater the chance of apical extrusion of the irrigant.^{67, 68} This must be avoided; were NaOCl to extrude past the apex, a catastrophic accident could occur.⁶⁹

_Manual-dynamic irrigation

Manual-dynamic irrigation involves gently moving a well-fitting gutta-percha master cone up and down in short 2 to 3 mm strokes within an instrumented canal, thereby producing a hydrodynamic effect and significant irrigant exchange.⁷⁰ Recent studies have shown that this irrigation technique is significantly more effective than automated-dynamic irrigation and static irrigation.^{9, 71, 72}

2

0

1

3



GREATER NEW YORK DENTAL MEETING

SCIENTIFIC MEETING:

Friday, November 29 -
Wednesday, December 4

EXHIBIT DATES:

Sunday, December 1 -
Wednesday, December 4

NO PRE-REGISTRATION FEE

Attend At No Cost

Never a pre-registration fee at the
Greater New York Dental Meeting

Mark Your Calendar

Educational Programs:
Friday through Wednesday,
November 29 - December 4

Exhibits:

Sunday through Wednesday,
December 1 - December 4

More than 600 Exhibitors

Jacob K. Javits Convention Center
11th Avenue between 34-39th Streets (Manhattan)

Headquarters Hotel

New York Marriott Marquis Hotel

Live Dentistry Arena - No Tuition

Latest Dental Technology & Scientific Advances

More Than 350 Scientific Programs

Seminars, Hands-on Workshops, Essays & Scientific
Poster Sessions as well as Specialty and Auxiliary
Programs

Educational Programs in various languages

Social Programs for the Entire Family

ENJOY NEW YORK CITY AT ITS BEST DURING THE
MOST FESTIVE TIME OF THE YEAR!

The Largest Dental Meeting/Exhibition/ Congress in the United States

89th ANNUAL SESSION




For More Information:

Greater New York Dental Meeting®

570 Seventh Avenue - Suite 800, New York, NY 10018 USA

Tel: (212) 398-6922 / Fax: (212) 398-6934

E-mail: victoria@gnydm.com / Website: www.gnydm.com

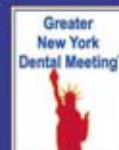

WWW.GNYDM.COM

Sponsored by: The New York County Dental Society and The Second District Dental Society

Machine-assisted agitation systems

Sonic irrigation

Sonic activation has been shown to be an effective method for disinfecting root canals, operating at frequencies of 1–6 kHz.^{73,74} There are several sonic irrigation devices on the market. The Vibringe allows delivery and sonic activation of the irrigating solution in one step. It employs a two-piece syringe with a rechargeable battery. The irrigant is sonically activated, as is the needle that attaches to the syringe. The EndoActivator is a more recently introduced sonically driven canal irrigation system.^{9,75} It consists of a portable handpiece and three types of disposable polymer tips of different sizes. The EndoActivator has been reported to effectively clean debris from lateral canals, remove the smear layer, and dislodge clumps of biofilm within the curved canals of molar teeth.⁹

Ultrasonics

Ultrasonic energy produces higher frequencies than sonic energy but low amplitudes, oscillating at frequencies of 25–30 kHz.^{9,76} Two types of ultrasonic irrigation are available. The first type is simultaneous ultrasonic instrumentation and irrigation, and the second type is referred to as passive ultrasonic irrigation operating without simultaneous irrigation (PUI). The literature indicates that it is more advantageous to apply ultrasonics after completion of canal preparation rather than as an alternative to conventional instrumentation.^{9,20,77} PUI irrigation allows energy to

be transmitted from an oscillating file or smooth wire to the irrigant in the root canal by means of ultrasonic waves.⁹ There is consensus that PUI is more effective than syringe needle irrigation at removing pulpal tissue remnants and dentine debris.^{78–80} This may be due to the much higher velocity and volume of irrigant flow that are created in the canal during ultrasonic irrigation.^{9,81} PUI has been shown to remove the smear layer; there is a large body of evidence with different concentrations of NaOCl.^{9,80–84} In addition, numerous investigations have demonstrated that the use of PUI after hand or rotary instrumentation results in a significant reduction in the number of bacteria,^{9,85–87} or achieves significantly better results than syringe needle irrigation.^{9,84,88,89}

Studies have demonstrated that effective delivery of irrigants to the apical third can be enhanced by using ultrasonic and sonic devices that demonstrate acoustic micro-streaming and cavitation.^{79,81,90,91} Acoustic micro-streaming is defined as the movement of fluids along cell membranes, which occurs as a result of the ultrasound energy creating mechanical pressure changes within the tissue. Cavitation is defined as the formation and collapse of gas and vapour-filled bubbles or cavities in a fluid.

The Apical Vapor Lock theory, proven *in vitro* by Tay, has been clinically demonstrated⁹² to also include the middle third by Vera: "The mixture of gases is originally trapped in the apical third, but then it might grow quickly by the nucleation of the smaller bubbles, forming a gas column that might not only impede penetration of the irrigant into the apical third but also push it coronally after it has been delivered into the canal." However, more recently Munoz⁹³ demonstrated that both: passive ultrasonic irrigation (PUI) and EndoVac are more effective than the conventional endodontic needle in delivering irrigant to WL of root canals."

This begs the efficacy question. Two recently published studies examined this issue with both systems by testing their ability to eliminate microorganisms during clinical treatment from infected root canal systems.^{94,95} Paiva found that after a supplementary irrigation procedure using PUI with NaOCl that 23% of the samples produced positive cultures. Cohenca's study examining the clinical efficacy of the EndoVac found no microbial growth either after post instrumentation irrigation or at the one⁹² week obturation appointment.

When questioning these diverse results one must remember that microbial hydrolysis via NaOCl is an equilibrium reaction. Hand demonstrated that a 50% reduction of NaOCl concentration resulted in a 300% reduction in dissolution activity. Accord-

Fig. 2_EndoVac set-up.

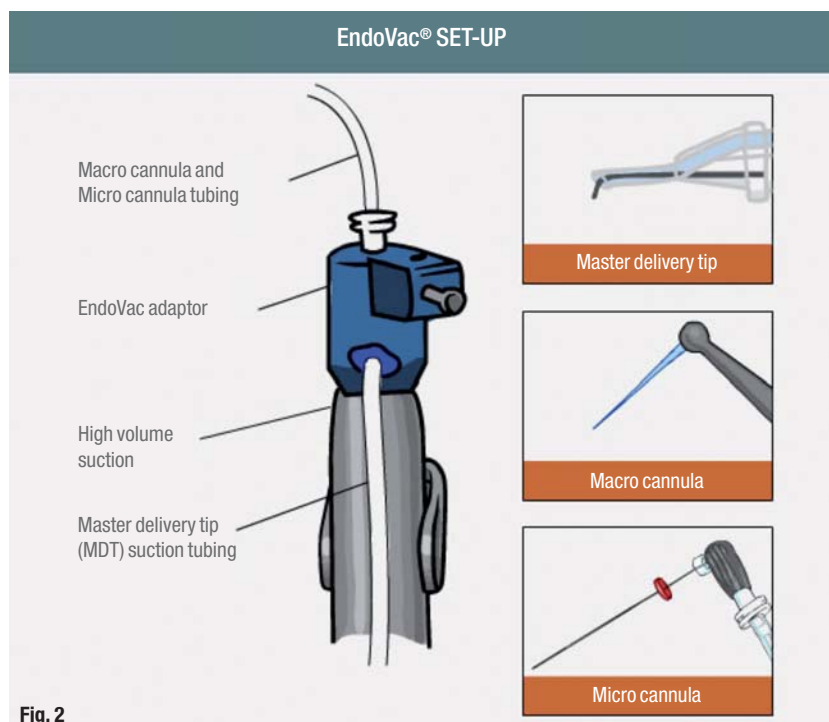


Fig. 2

ingly, one must consider both the delivery of the irrigant to full working length, via PUI or apical negative pressure and the total volume of NaOCl exchanged. The volume of an instrumented root canal 19 mm long shaped to a #35 with a 6% instrument equals .014 cc. Paiva described placement of NaOCl via a NaviTip (ULTRADENT) at WL – 4 mm during instrumentation and discussed using PUI with #15 K-file at WL – 1 mm. Prior to PUI, 2 ml of NaOCl was injected into the canal; however, this could not have filled the apical four millimeters⁹⁵ due to the apical vapor lock. According to Munoz, the canal was most likely immediately filled with ultrasonically activated NaOCl for one minute⁹², but as just described – only about .014cc would have been effectively available for this exchange and activation. In contrast, the Apical Negative Pressure protocol described by Cohenca et al. approximately 2 ml of NaOCl actively passes through the complete WL for one⁹² minute.⁹⁶ The difference in volumetric exchange equals $2/.014 = 14,200\%$ and likely explains the disinfection differential.

_The plastic rotary F File

Although sonic or ultrasonic instrumentation is more effective at removing residual canal debris than rotary endodontic files are,¹⁰⁴ and irrigation solutions are often unable to remove this during endodontic treatment, many clinicians still do not incorporate it into their endodontic instrument armamentarium. The common reasons given for not using sonic or ultrasonic filing are that it can be time-consuming to set up, an unwillingness to incur the cost of the equipment, and lack of awareness of the benefits of this final instrumentation step in endodontic treatment.

It is for these reasons that an endodontic polymer-based rotary finishing file was developed. This new, single-use, plastic rotary file has a unique file design with a diamond abrasive embedded into a non-toxic polymer. The F File will remove dentinal wall debris and agitate the NaOCl without enlarging the canal further.

_Pressure-alternation devices

Rinsendo irrigates the canal by using pressure-suction technology. Its components are a handpiece, a cannula with a 7 mm exit aperture, and a syringe carrying irrigant. The handpiece is powered by a dental air compressor and has an irrigation speed of 6.2 ml/min. Research has shown that it has promising results in cleaning the root-canal system, but more research is required to provide scientific evidence of its efficacy. Periapical extrusion of irrigant has been reported with this device.^{101, 102}

_The EndoVac apical negative-pressure system

The EndoVac apical negative-pressure irrigation system has three components: the Master Delivery Tip, MacroCannula and MicroCannula. The Master Delivery Tip simultaneously delivers and evacuates the irrigant (Fig. 2). The MacroCannula is used to suction irrigant from the chamber to the coronal and middle segments of the canal. The MacroCannula or MicroCannula is connected via tubing to the high-speed suction of a dental unit. The Master Delivery Tip is connected to a syringe of irrigant and the evacuation hood is connected via tubing to the high-speed suction of a dental unit.⁵⁶ The plastic MacroCannula has an open end of ISO size 0.55 mm in diameter with a 0.02 taper and is attached to a handpiece for gross, initial flushing of the coronal and mid-length parts of the root canal. The MicroCannula contains 12 microscopic holes and is capable of evacuating debris to full working length.¹⁰² The ISO size 0.32 mm diameter stainless-steel MicroCannula has four sets of three laser-cut, laterally positioned offset holes adjacent to its closed end, 100 μ in diameter and spaced 100 μ apart. This is attached to a finger piece for irrigation of the apical part of the canal when it is positioned at working length. The MicroCannula can be used in canals that are enlarged with endodontic files to ISO size 35.04 or larger.

During irrigation, the Master Delivery Tip delivers irrigant to the pulp chamber and siphons off the excess irrigant to prevent overflow. Both the MacroCannula and MicroCannula exert negative pressure that pulls fresh irrigant from the chamber, down the canal to the tip of the cannula, into the cannula, and out through the suction hose. Thus, a constant flow of fresh irrigant is delivered by negative pressure to working length. A recent study showed that the volume of irrigant delivered was significantly higher than the volume delivered by conventional syringe needle irrigation within the same period,⁴⁶ and resulted in significantly more debris removal at 1 mm from working length than did needle irrigation. During conventional root-canal irrigation, clinicians must be careful when determining how far an irrigation needle is placed into the canal. Recommendations for avoiding NaOCl incidents include not binding the needle in the canal, not placing the needle close to working length, and using a gentle flow rate when using positive-pressure irrigation.¹⁰³ With the EndoVac, in contrast, irrigant is pulled into the canal at working length and removed by negative pressure. Apical negative pressure has been shown to enable irrigants to reach the apical third and help overcome apical vapour lock.^{46, 104} In addition, with respect to isthmus cleaning, although it is not possible to reach and clean the isthmus area with instruments, it is not impossible to reach and thoroughly clean these



Fig. 3

Fig. 3 Irrigation accident with widespread trauma.

areas with NaOCl when the method of irrigation is safe and efficacious. In studies comparing the EndoActivator,¹⁰⁵ passive ultrasonic,¹⁰⁵ the F File,¹⁰⁵ the manual-dynamic Max-i-Probe (DENTSPLY Rinn),^{105, 106} the Pressure Ultrasonic¹¹¹ and the EndoVac,¹⁰⁶ only the EndoVac was capable of cleaning 100% of the isthmus area.

Apart from being able to avoid air entrapment, the EndoVac system is also advantageous in its ability to deliver irrigants safely to working length without causing their undue extrusion into the periapex,^{46, 102} thereby avoiding NaOCl incidents. It is important to note that it is possible to create positive pressure in the pulp canal if the Master Delivery Tip is misused, which would create the risk of a NaOCl incident. The manufacturer's instructions must be followed for correct use of the Master Delivery Tip.

Sodium hypochlorite incidents

Although a devastating endodontic NaOCl incident is rare,¹⁰⁷ the cytotoxic effects of NaOCl on vital tissue are well established.¹⁰⁸ The associated seque-

lae of NaOCl extrusion have been reported to include life-threatening airway obstructions,¹⁰⁹ facial disfigurement requiring multiple corrective surgical procedures,¹¹⁰ permanent paraesthesia with loss of facial muscle control,⁶⁹ and—the least significant consequence—tooth loss.¹¹¹

Although the exact aetiology of the NaOCl incident is still uncertain, based on the evidence from actual incidents and the location of the associated tissue trauma, it would appear that an intravenous injection may be the cause. The patient shown in Figure 3 demonstrates a widespread area of tissue trauma that is in contrast to the characteristics of NaOCl incident trauma reported by Pashley.^{108, 112} This extensive trauma, and particularly involving the pattern of ecchymosis around the eye, could only have occurred if the NaOCl had been introduced intravenously to a vein close to the root apex through which extrusion of the irrigant occurred and the irrigant then found its way into the venous complex. This would require positive pressure apically that exceeded venous pressure (10 mg of Hg). In one *in vitro* study, which used a positive-pressure needle irrigation technique to mimic clinical conditions and techniques, the apical pressure generated was found to be eight times higher than the normal venous pressure.¹¹³

This does not imply that NaOCl can or should be excluded as an endodontic irrigant; in fact, its use is critical, as has been discussed in this article. What this does imply is that it must be delivered safely.

Safety first

In order to compare the safety of six current intra-canal irrigation delivery devices, an *in vitro* test was conducted using the worst-case scenario of apical extrusion, with neutral atmospheric pressure and an open apex.¹⁰² The study concluded that the EndoVac did not extrude irrigant after deep intra-canal delivery and suctioning of the irrigant from the chamber to full working length, whereas other devices did. The EndoActivator extruded only a very small volume of irrigant, the clinical significance of which is not known.

Mitchell and Baumgartner tested irrigant (NaOCl) extrusion from a root canal sealed with a permeable agarose gel.¹¹⁴ Significantly less extrusion occurred using the EndoVac system compared with positive-pressure needle irrigation. A well-controlled study by Gondim et al. found that patients experienced less post-operative pain, measured objectively and subjectively, when apical negative-pressure irrigation was performed (EndoVac) than with apical positive-pressure irrigation.¹¹⁵

_Efficacy

In vitro and *in vivo* studies have demonstrated greater removal of debris from the apical walls and a statistically cleaner result using apical negative-pressure irrigation in closed root-canal systems with sealed apices. In an *in vivo* study of 22 teeth by Siu and Baumgartner, less debris remained at 1 mm from working length using apical negative pressure compared to use of traditional needle irrigation, while Shin et al. found in an *in vitro* study of 69 teeth comparing traditional needle irrigation with apical negative pressure that these methods both resulted in clean root canals, but that apical negative pressure resulted in less debris remaining at 1.5 and 3.5 mm from working length.^{46, 104, 116} When comparing root-canal debridement using manual-dynamic agitation or the EndoVac for final irrigation in a closed system and an open system, it was found that the presence of a sealed apical foramen adversely affected debridement efficacy when manual-dynamic agitation was used, but did not adversely affect results when the EndoVac was used. Apical negative-pressure irrigation is an effective method to overcome the fluid-dynamic challenges inherent in closed root-canal systems.¹¹⁷

_Microbial control

Hockett et al. tested the ability of apical negative pressure to remove a thick biofilm of *E. Faecalis*, finding that these specimens rendered negative cultures obtained within 48 hours, while those irrigated using traditional positive-pressure irrigation were positive at 48 hours.⁹⁹

One study found that apical negative-pressure irrigation resulted in similar bacterial reduction to use of apical positive-pressure irrigation and a triple antibiotic in immature teeth.¹¹⁸ In a study comparing the use of apical positive-pressure irrigation and a triple antibiotic that has been utilised for pulpal regeneration/revascularisation in teeth with incompletely formed apices (Trimix = Cipro, Minocin, Flagyl) versus use of apical negative-pressure irrigation with NaOCl, it was found that the results were statistically equivalent for mineralised tissue formation and the repair process.¹¹⁹ Using apical negative pressure and NaOCl also avoids the risk of drug resistance, tooth discoloration, and allergic reactions.^{120, 121}

_Conclusion

Since the dawn of contemporary endodontics, dentists have been syringing NaOCl into the root-canal space and then proceeding to place endodontic instruments down the canal in the belief that

they were carrying the irrigant to the apical termination. Biological, scanning electron microscopy, light microscopy, and other studies have proven this belief to be in error. NaOCl reacts with organic material in the root canal and quickly forms micro-bubbles at the apical termination that coalesce into a single large apical vapour bubble with subsequent instrumentation. Since the apical vapour lock cannot be displaced via mechanical means, it prevents further NaOCl flow into the apical area. The safest method yet discovered to provide fresh NaOCl safely to the apical terminus to eliminate the apical vapour lock is to evacuate it via apical negative pressure. This method has also been proven to be safe because it always draws irrigants to the source via suction—down the canal and simultaneously away from the apical tissue in abundant quantities.¹²² When the proper irrigating agents are delivered safely to the full extent of the root-canal terminus, thereby removing 100% of organic tissue and 100% of the microbial contaminants, success in endodontic treatment may be taken to levels never seen before._

Editorial note: A complete list of references is available from the publisher.

This article has been reprinted in part from G. Glassman, Safety and Efficacy Considerations in Endodontic Irrigation (PenWell, January 2011).

_about the author

roots



Dr Gary Glassman graduated from the University of Toronto Faculty of Dentistry in 1984 and graduated from the Endodontology Program at Temple University in 1987, where he received the Louis I. Grossman Study Club Award for

academic and clinical proficiency in endodontics. The author of numerous publications, he lectures globally on endodontics and is on the staff at the University of Toronto Faculty of Dentistry in the Graduate Department of Endodontics. He is a fellow of the Royal College of Dentists of Canada, and the endodontic editor for the Oral Health journal. He maintains a private practice, Endodontic Specialists, in Toronto, Ontario, Canada. He can be reached through his website, www.rootcanals.ca