

Er:YAG laser scanner for **implant site** preparation

Drs Norberto Berna, Giovanni Olivi, Luca Marigo & Massimo Cordaro, Italy

Introduction

Numerous studies have reported the ability of mid-infrared lasers (Er:YAG laser with a wavelength of 2,940nm and Er,Cr:YSGG laser with a wavelength of 2,780nm) to ablate hard biological tissue without creating thermal damage.¹⁻⁵ Among the various applications, the possibility of preparing the implant site entirely with laser irradiation was originally illustrated *in vivo* by Dr Berna, one of the authors of this article, in 2003. A patent for this form and method was issued and registered at that time. The first group of 62 patients received this method of treatment, using an Er,Cr:YSGG laser, between November 2001 and December 2002. Since then, another larger group of patients have been treated with this method. An Er:YAG laser and a dedicated laser scanner handpiece have never been used in human patients for this proce-

cedure, however. The major obstacles to the wide use of this technique are the time to create the osteotomy,⁶ the low energy capability of the devices previously available and the time required to learn to use a handpiece that works in no contact. In the present study, a high performance Er:YAG laser was used with a tipless laser scanner handpiece that allows more precision, a higher energy output and a shorter pulse duration in comparison with the previously used device.

Materials and methods

A 67-year-old patient, a non-smoker without any systemic diseases, was examined using CBCT (PaX-i3D Smart, Vatech) in order to assess the surgical area, the bone volume and the bone density in the edentulous area, region #15 (Fig. 1). Before the surgery, the patient

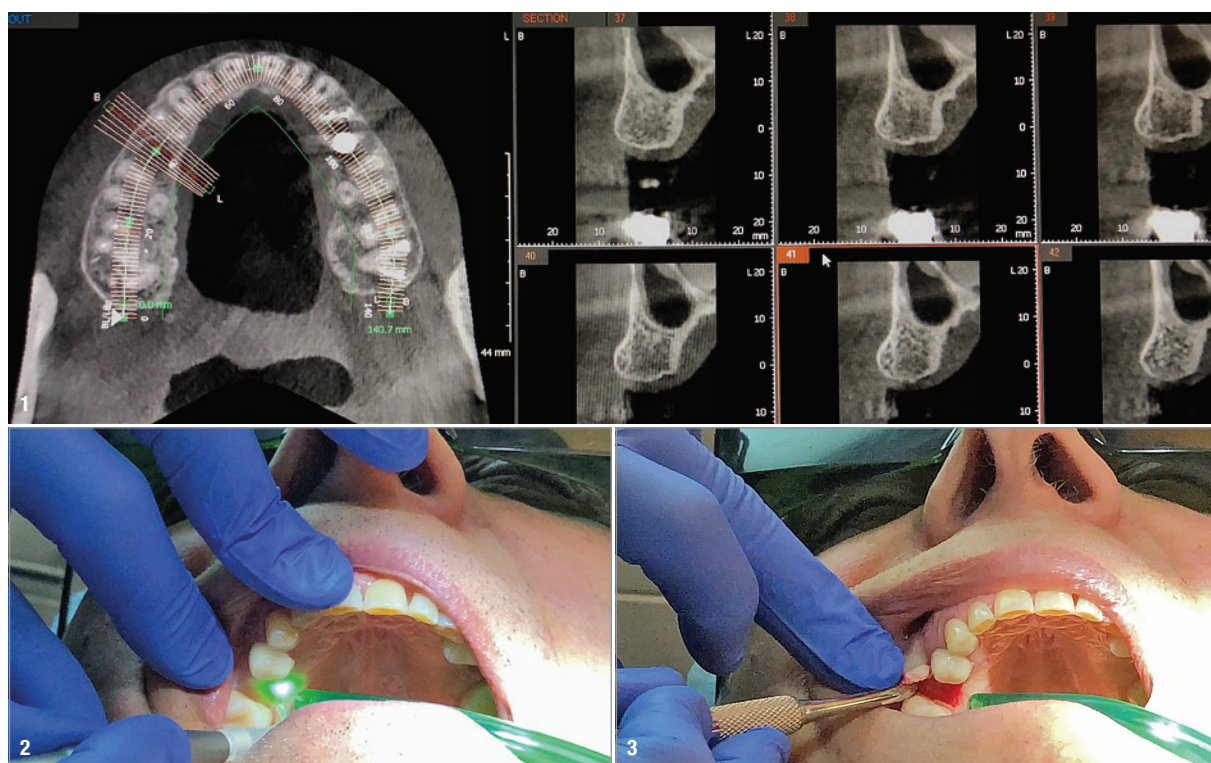


Fig. 1: A CBCT scan was taken pre-op. **Fig. 2:** Flap incision using a LightWalker Nd:YAG laser head (200 µm diameter fibre, 3W, 70Hz, MSP). **Fig. 3:** Flap reflection was carried out with a Prichard elevator.

EXCEPTIONALLY NO ADVERTISING AT THIS POINT, BUT A GREAT THANK YOU TO EVERYONE WHO HAS RESEARCHED AND DEVELOPED WITH US IN THE PAST YEARS.

Henning Alberti, h4sterne.de, Visual Designer and Photographer | **Jandar Amir**, Dentallabor Moss GmbH, Trainee | Prof. Dr. med. habil. Dipl.-Ing. **Rainer Bader**, Management, Orthopädische Klinik und Poliklinik Universitätsmedizin Rostock, Forschungslabor für Biomechanik und Implantattechnologie FORBIOMIT | Dr. **Angela Bergmann**, pipVerlag, Medical Writing & Editorial Staff | Dr. med. Dr. med. dent. **Wilken Bergmann**, Bundeswehrkrankenhaus Hamburg, Specialist in oral and maxillofacial surgery Focus on implantology | Dr. med. dent. **Michael Berthold Seehofer**, Zahnärzte in München, Oral Surgery Specialist | **Nico Bizer**, NB VIDEOFACTORY, Managing Director | Family **Bjelajac** | **Ralf Blank**, TBI Technologie-Beratungs-Institut GmbH, Managing Director | Dr. **Catrin Bludszuweit-Philipp**, ASD Advanced Simulation & Design GmbH, Managing Director | **Uwe Bräuer**, GENIUS Venture Capital GmbH, Managing Director | **Ulrike Burmeister**, Poliklinik für Zahnerhaltung und Parodontologie, Universitätsmedizin Rostock, Senior Consultant | Dr. med. dent. **Claudio Cacaci**, IMPLANTAT COMPETENCE CENTRUM, Gemeinschaftspraxis Cacaci & Randelzhofer, Dentist & Oral Surgeon | **Catrin Crasemann**, Creative Consultant | Dr. med. Dr. med. dent. **Michael Dau**, Klinik und Poliklinik für Mund-, Kiefer und Plastische Gesichtschirurgie, Specialist in Oral and Maxillofacial Surgery, Universitätsmedizin Rostock | **Denise Nielsen**, TIZIO Hybrid Implants GmbH, Sales & Logistic | **Mario Dombrowski**, FMZ GmbH, Head of Manufacturing | PD. Dr. rer. nat. **Nadja Engel**, Klinik und Poliklinik für Mund-, Kiefer und Plastische Gesichtschirurgie, Head of Research Laboratory | **Stefanie Flentje**, Dentallabor Moss GmbH, Master of Dental Technician | **Andre Flentje**, grafjik Kommunikationsdesign, Graphic Design and Photographer | **Karin Franke**, Schackow Joksch Rechtsanwälte PartG mbB, Attorney | Prof. Dr. med. Dr. med. dent. **Bernhard Frerich**, Klinik und Poliklinik für Mund-, Kiefer und Plastische Gesichtschirurgie, Universitätsmedizin Rostock, Clinic Director | **Rene Friedrich**, ZM Präzisionsdentaltechnik GmbH, Dental Technican | Dr. **Andreas Frost**, Dental Office, Dentist | Dr. **Claudia Fruijtter-Pölloth**, CATS Consultants GmbH, Managing Director | Dr. **Antonia Fuchs**, Dental Office, Dentist | Dr. **Cornelia Ganz**, TIZIO Hybrid Implants GmbH, Head of Research & Development | **Reinhard Gedack**, D.I.E. WERBEFILM GmbH, Senior Executive Producer | **Patrik Gölz**, paddefilm, Best Cameraman | Dr. **Gian-Garlo Gullo**, Metoxit AG, Managing Director | Dipl.-Ing. **Dierk Hansmann**, Patentanwälte Hansmann & Hansmann, Patent Attorney | **Robert Hauschild**, FMZ GmbH, Research & Development | Dr. med. Dr. med. dent. **Anders Henningsen**, ELBE MKG, Specialist in Oral and Maxillofacial Surgery Focus on Implantology | **Ulrike Hinz**, TIZIO Hybrid Implants GmbH, Quality Management Officer | Dr. med. **Michael Hopp**, Dental Office, Dentist | Dr. med. dent. **Ulrike Hoppe**, Dentalmedizin Hamburg, Dentist | **Frank Idel**, Gebr. Brasseler GmbH & Co. KG, Senior Manager Engineering & Quality Control | **Hubert Joksche**, Schackow Joksch Rechtsanwälte PartG mbB, Attorney | Dr. rer. hum. Dipl.-Biol. **Anika Jonitz-Heincke**, Forschungslabor für Biomechanik und Implantattechnologie FORBIOMIT, Orthopädische Klinik und Poliklinik Universitätsmedizin Rostock, Head of the Tissue Regeneration Group | **Stephan Jung**, Dentallabor Moss GmbH, Dental Technican | **Kerstin Jung**, Kommunikation Dental, Editor | **Rita Kohlbach**, FMZ GmbH, Logistic | Prof. Dr. med. Dr. rer. nat. **Markus Kipp**, Institut für Anatomie, Universitätsmedizin Rostock, Institute Director | **Wolfgang Klinger**, ACF Amberger Central Fräs-Center GmbH, Managing Director | Priv.-Doz. Dr.-Ing. habil. **Daniel Klüß**, innoproof GmbH, Managing Director | **Thomas Knüppel**, KLOPSCH & PARTNER RECHTSANWÄLTE mbB, Attorney | **Susanne Knüppel**, STEUERBOARD, Tax Advice & Economic Advice | Dr. **Stefan Köbel**, Decema GmbH, Managing Director | Dipl.-Ing. **Uwe Koch**, primec GmbH, Managing Director | **Joanna Moss**, Dentallabor Moss GmbH, Trainee | Prof. **Hans-Georg Neumann**, DOT GmbH, Managing Director | Dr. **Barbara Neumann**, Composita GmbH | **Hans-Christian Neumann**, Composita GmbH, Managing Director | Dr. med. dent. **Thomas Niedermeier**, Dental Office, Dentist & Implantologist | **Cord Obenhausen**, Gebr. Brasseler GmbH & Co. KG, Regional Sales Manager | Dr. **Hartmut Ohm**, ABZ Abrechnungs- und Beratungsgesellschaft für Zahnärzte eG, Chairman of the Board | **Iris Peters**, STEUERBOARD, Tax Advice & Economic Advice | Dr. **Pantelis Petrakakis**, pipVerlag, Freelance Journalist & Dentist | Dr. **Andreas Pippig**, Dental Office, Dentist | **Otto Prandtner**, Plattform für feinste Dentaltechnologie GbR, Dental Technican | Dr. med. dent. **Friedrich Preusse**, Dres. Kleier & Schlieper, Oral Surgery Specialist | **Dirk Pfützner**, FMZ GmbH, Managing Director | **Katja Quednau**, FMZ GmbH, Quality Management | **Lars Quiring**, D.I.E. WERBEFILM GmbH, Compositing Artist | **Ronny Raasch**, OstseePharma Consulting GmbH, Managing Director | Dr. med. dent. **Peter Randelzhofer**, IMPLANTAT COMPETENCE CENTRUM, Gemeinschaftspraxis Cacaci & Randelzhofer, Dentist | **Daniel Reimer**, FMZ GmbH, Production | Dr. med. dent. **Morten Reimer**, Zahnarztpraxis, Dentist | Dr. **Michael Reise**, estetic ceram ag, Head of Research & Development | Dr. **Bärbel Riemer-Krammer**, Dental Office, Dentist | **Nicole Sabielny**, TIZIO Hybrid Implants GmbH, Quality Management Officer | Dr. med. dent. **Reza Saeidi-Pour**, Seehofer – Zahnärzte in München, Dentist | **Nico Schade**, STEUERBOARD, Tax Advice & Economic Advice | **Julian Scheinkönig**, SCHEINKÖNIG, Creativ Director | **Carsten Schiele**, Rostock Business and Technology Development GmbH, Project Manager | Dr. **Benedikt Seeber**, Metoxit AG, CTO | Dr. med. dent. **Peter Seehofer**, Seehofer – Zahnärzte in München, Dentist | **Katrin Seidel**, TIZIO Hybrid Implants GmbH, Office | Prof. Dr.-Ing. **Hermann Seitz**, Lehrstuhl für Mikrofluidik, Universität Rostock, Chair Director | **Marc Simon**, D.I.E. WERBEFILM GmbH, Executive Producer | **Lara Speitmann**, projekt rk Messe Design Werbung, Project Manager | **Egbert Sprenger**, estetic ceram ag, Managing Director | Prof. Dr. Dr. h.c. **Frank Stein**, Healthcare Projects Consulting & Management, QM Consulting | **Marianne Steinbeck**, pipVerlag, Publisher | Dr. **Bernhard Steinfeldt**, Zahnarztpraxis, Dentist | **Anja Stolzenburg**, FMZ GmbH, Logistic | Prof. Dr. **Joachim Tinschert**, Universitätsklinikum Aachen, Klinik für Zahnärztliche Prothetik, Dentist | **Nico Tjarks**, Homepage Helden GmbH, Managing Director | **Thomas Ulrich**, ulrich-fotodesign, Photographer | **Danny Vogel**, Forschungslabor für Biomechanik und Implantattechnologie FORBIOMIT, Orthopädische Klinik und Poliklinik Universitätsmedizin Rostock, Head of the Experimental Biomechanics Group | **Martin Vollbrecht**, Dental Balance GmbH, Managing Director | **Helge Vollbrecht**, Dental Balance GmbH, Managing Director | Prof. Dr. med. habil. **Brigitte Vollmar**, Universitätsmedizin Rostock, Institute Director Experimental Surgery | Dr. **Wolfram Weber**, Decema GmbH, Managing Director | **Christian Weiß**, Rostock Business and Technology Development GmbH, Managing Director | **Wolfgang Weisser**, CTB-Zirkel, Master of Dental Technician | **Christina Wels**, ZM Präzisionsdentaltechnik GmbH, Dental Technican | **Steve Wendland**, Rostock Business and Technology Development GmbH, Senior Consultant Technology Advancement | **Daniel Wolter**, Klinik und Poliklinik für Mund-, Kiefer und Plastische Gesichtschirurgie, Medical-Technical Assistance | **Kastriot Xhoxhaj**, estetic ceram ag, Sales and Marketing | **Ralf Zerbe**, FMZ GmbH, Production | Dr. rer. hum. Dipl.-Ing. **Carmen Zietz**, innoproof GmbH, Managing Director | Dr. med. dent. **Catarina Zimmermann**, Dentalmedizin Hamburg, Dentist & Oral Surgeon | **Elbpäntent Marschall & Partner** Rechtsanwaltskanzlei, Patent Attorney | **Claudia Voss**, TIZIO Hybrid Implants GmbH, Sales Director | **Jennifer Wilken**, TIZIO Hybrid Implants GmbH, Managing Director



TIZIO H6

TIZIO F3

WE WOULD BE HAPPY TO PRESENT OUR INNOVATIVE
HYBRID TECHNOLOGY TO YOU IN DETAIL. YOU CAN
FIND CONTACT AND FURTHER INFORMATION AT
WWW.TIZIOIMPLANTS.COM
WE ARE PLEASED TO MEET YOU. STAY HEALTHY.



TIZIO HYBRID IMPLANTS
The Root of the Future



Fig. 4: The stabilised laser preparing the osteotomy, operated at 380mJ, 20Hz and 50 microseconds. **Fig. 5:** The Osstell handpiece detecting primary stability. **Fig. 6:** Determining the ISQ value. **Fig. 7:** Clinical situation post-op. **Fig. 8:** Radiograph taken immediately after surgery. **Fig. 9:** Radiograph taken after eight months of loading.

had received all the information regarding the treatment and the possible alternative treatment through a personalised informed consent form. The implant insertion axis has been planned for the best functional result of the prosthesis. A customised dental resin holder was created to support the laser scanner handpiece intra-orally, in the correct position according to the insertion axis of the implant. Local anaesthesia was administered using articaine (1 : 100,000). A full-thickness incision was performed using an Nd:YAG laser (1,064nm wavelength; 200µm diameter fibre; MSP: 3W, 70Hz; LightWalker AT, Fotona) on the palatal paramedian line; two mesial and distal releasing incisions were also performed without involving the papillae (Fig. 2). The access flap was then reflected with a Prichard elevator (Fig. 3).

An Er:YAG laser with a wavelength of 2,940nm (LightWalker AT) equipped with a laser scanner handpiece (X-Runner, Fotona) was used (Fig. 4). The laser parameters used were 380mJ and 20Hz, delivered with a super-short pulse (50 microseconds). An external source of sterile saline solution stored at 5°C in a refrigerator was used, and the saline was delivered via a peristaltic pump to promote photothermal ablation and to reduce the temperature in the surgical site. The scanner allows one to program and precisely perform a circular osteotomy of 3.5mm in diameter, the same diameter as the implant manufacturer's final drill. During the osteotomy, the insertion depth was checked using a millimetre probe, until the preset depth of 12mm was reached. The author prefers to place implants 2mm deeper sub-crestally to prevent angular resorption and to manage the emergence profile of the prosthesis more effectively. A tapered screw implant made of Grade IV titanium and with a sandblasted and acid-etched surface (HELI, IDC) was inserted. It had a maximum diameter of 4.2mm on the external thread and a length of 10.0 mm. Once inserted, the implant stability quotient (ISQ) was determined using the Osstell handpiece (Osstell; Figs. 5 & 6). The flap was sutured (Fig. 7), a postoperative radiograph was taken (Fig. 8), and after five days, the sutures were removed. At that time, the patient was asked to assess the postoperative pain he had experienced by assigning a numeric value of between 0 and 10 on a verbal numeric scale.



Results

The total clinical time for preparing the osteotomy was approximately 7 minutes. The implant had a high primary stability value at the time of insertion: the ISQ score measured in the buccolingual direction was 84 and the score measured in the mesiodistal direction was 81. The reported numeric value of postoperative pain was 1. At the second stage of the implant treatment, which was performed after 40 days, new ISQ values of 84 buccolingually and 82 mesiodistally were determined. After three months, the values had increased to 86 buccolingually

SkyPulse®

A new generation of compact
Fotona lasers

A compact, easy
to use Er:YAG Laser
tailored for Your
practice!

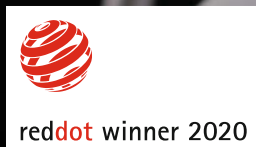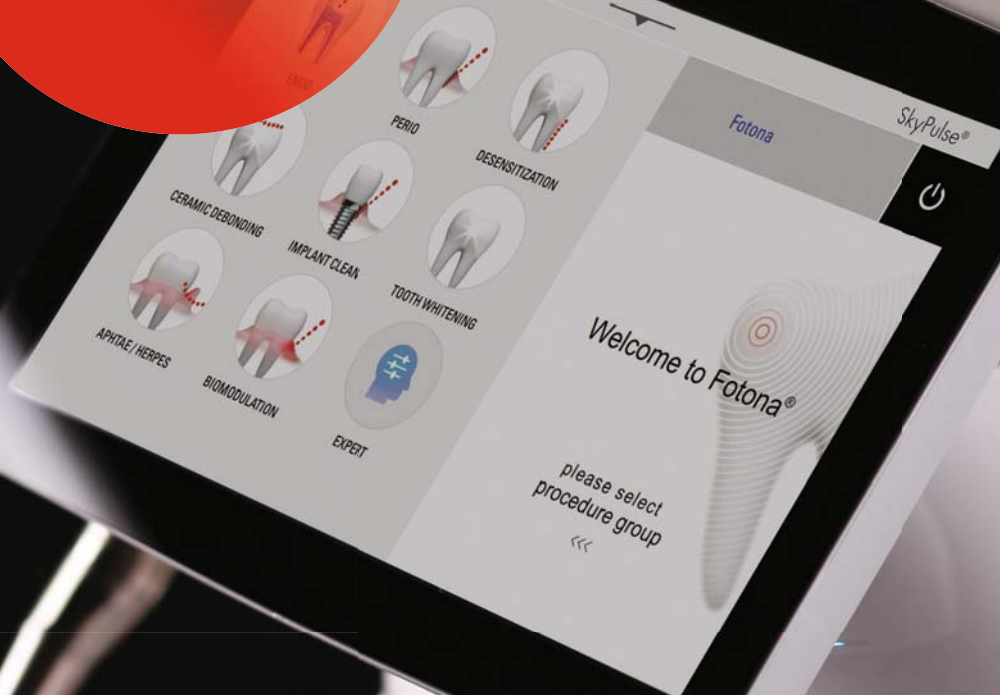

Committed to Engineering

The Highest Performance, Best Made Laser Systems in the World

Fotona®
choose perfection

www.fotona.com

and 84 mesiodistally. The radiographic control image taken after eight months of loading showed complete osseointegration of the implant with very good healing of the bone around the neck of the abutment (Fig. 9).

Discussion

The Er:YAG laser technology used in this study allowed a reduction of bone preparation time compared with the laser system used previously. Also considering the absence of pauses needed to replace the drills, the total time for preparing the osteotomy was very close to that of the conventional method.⁷ However, there are other advantages of using a laser for preparing osteotomies. Besides a comparable preparation time, the higher energy delivered through a tipless stabilised laser scanner handpiece to the target, with a super-short pulse and efficient water spray, reduced the minor thermal damage to the bone cells adjacent to the osteotomy site and avoided the smear layer typically associated with conventional drill preparation.^{8,9} This could reduce the time necessary for osseointegration.¹⁰ The laser site preparation is very precise and compatible with different implant geometries, owing to the absence of micro-movements and skidding induced by the impact of drills on the bone, especially for large-diameter final drills.¹¹ Such a scanner technology is of major significance for developing laser use in implantology further, as it allows osteotomies of different sizes and shapes (circular and possibly oval) to be created with very high precision and compatibility with new implant geometries.

The ISQ score analysis used in this clinical study reports the values of resonance frequency analysis using the Osstell handpiece. The scale of values ranges from 0 to 100, and the lowest threshold value for single implant loading, even in cases of immediate placement, corresponds to 70. This is the only non-invasive method accepted by the international scientific community for clinical evaluation of osseointegration.¹²⁻¹⁵ In the clinical case described in this article, the high primary stability value increased postoperatively over the course of 40 days, indicating the absence of bone damage during the laser procedure and that cell healing had started promptly and successfully. Indeed, erbium laser bone irradiation has been reported to increase release of bone morphogenetic protein, suggesting quicker osseointegration and stability as a result.¹⁶ Among the other advantages of laser osteotomy preparation, the lack of contact, vibration and pressure made for a more comfortable experience for the patient during the intervention. Finally, the immediate postoperative period, as reported by the patient, was free of pain and complications.

Conclusion

This first clinical case highlights the many advantages of the use of the Er:YAG laser to perform an implant oste-

otomy without drills. The absence of postoperative pain and complications and the very high value of primary implant stability described in this case report suggest that the use of an Er:YAG laser and a dedicated laser scanner handpiece is a viable treatment alternative for creating osteotomies. These positive results support the realisation of a pilot study with a greater number of implants.

about the authors




Dr Norberto Berna graduated in Medicine and Surgery and specialised in dentistry and stomatology. He holds an MSc in Odontology and Forensic Science. Since 1987, he has been working in his private office with lasers of different wavelengths. He holds three patents and is frequently invited to speak at international conferences.

He is a joint scientific coordinator of the master's degree in laser dentistry at the Università Cattolica del Sacro Cuore in Rome in Italy. Between 2012 and 2018, he was an adjunct professor at the University of Genoa in Italy.



Dr Giovanni Olivi graduated cum laude in Medicine and Surgery (MD) and in Dentistry (DDS) from universities in Rome. He obtained a postgraduate diploma in laser dentistry from the University of Florence in Italy in 2002 and master status from the Academy of Laser Dentistry in Coral Springs in Florida in the US in 2009. Since 2007,

Dr Olivi has been a contracted professor and a joint scientific coordinator of the master's degree in laser dentistry at the Università Cattolica del Sacro Cuore. He is the author of over 70 peer-reviewed articles, several chapters in dentistry textbooks and four books on laser dentistry.

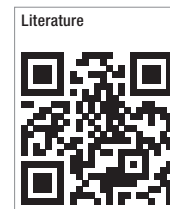

contact

Dr Norberto Berna

info@laserimplantology.com

Dr Giovanni Olivi

olivi@inlaser.it

Conical? Parallel?

copaSKY!

The unique hybrid connection
for sophisticated restorations!



Leading in Immediate Restoration – powered by physiological Prosthetic!

DENTAL INNOVATIONS
SINCE 1974

bre^{dent} group

bre^{dent} UK Ltd · Unit 26, Tapton Park Innovation Centre · Brimington Road · Chesterfield S41 0TZ · T: +44 1246 559599 · F: +44 1246 557777 · www.bre^{dent}.co.uk · @: enquiries@bre^{dent}.co.uk