


# Laser-enhanced endodontic treatment

**Author:** Dr Gregori M. Kurtzman, USA

**Endodontic success** is predicated on the ability to debride and clean the canal system. That canal system is a complex array of accessory and lateral canals, fins and other anatomical areas inaccessible to endodontic files (Fig. 1). As practitioners, we are able to clean the main canals with files, either hand or rotary. But they cannot mechanically remove pulpal tissue and debris from the canal anatomy present adjacent to the main canals.

Treatment success requires elimination of the pulpal tissue and associated bacteria from this anatomy, so that it can be sealed during the obturation phase of treatment. As only one thing can occupy a space at a time, obturation material cannot fill areas still occupied by pulpal tissue. Success is dependent on disinfection and debridement of the canal system so that it may be sealed during obturation. Irrigation has long been accepted as a key factor of treatment to achieve those goals.

**Fig. 1:** Anatomy of the canal system demonstrating accessory canals, fins and lateral canals which are not accessible with endodontic files as shown in cleared teeth.

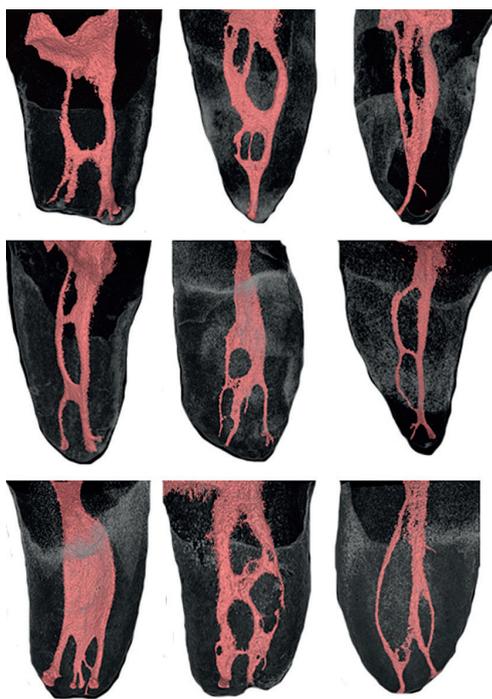


Yet, complete clearing of residual bacteria especially in the apical portion of the canal system has been difficult to achieve with traditional methods using even sodium hypochlorite (NaOCl) solutions (Fig. 2). Studies have demonstrated that the additional use of an Er:YAG laser to activate the irrigation solution greatly enhances not only the efficiency of the irrigation solutions advocated (NaOCl and EDTA) but also improves disinfection of the canal system, clearing accessory so that it may be sealed during obturation (Figs. 3 & 4).

## Irrigation the key to endodontic success

Although, instrumentation with files is important to enlarging the canals and ready them to be obturated, debris consisting of pulpal tissue and associated bacteria is not effectively removed by files. Irrigation with an appropriate solution is required to remove that debris from the canal walls. NaOCl is still the accepted irrigant due to its tissue dissolving ability and antibacterial nature. Yet, it cannot effectively reach far beyond the main canals to remove the residual tissue. Tissue dissolution can be enhanced to more effectively remove pulpal tissue/bacteria and also reach further into the accessory anatomy to allow better sealing of the canal system improving treatment success.

Smear layer within the canal system plays a factor in success in endodontic treatment. The smear layer contains bacteria which when left within the canal anatomy may lead to recurrence of infection endodontically. When compared to traditional irrigation methods, laser enhanced irrigation has demonstrated better intracanal smear layer removal.<sup>1</sup> As *Enterococcus faecalis* has been routinely linked to endodontic failures, and is a common occupant of the oral cavity, elimination of this bacteria is critical to prevention of reinfection of the canal system. NaOCl as an irrigant has not shown to be effective in elimination of *E. faecalis*, yet when combined with laser enhanced irrigation with NaOCl this bacteria has been eliminated in the canal anatomy.<sup>2</sup>

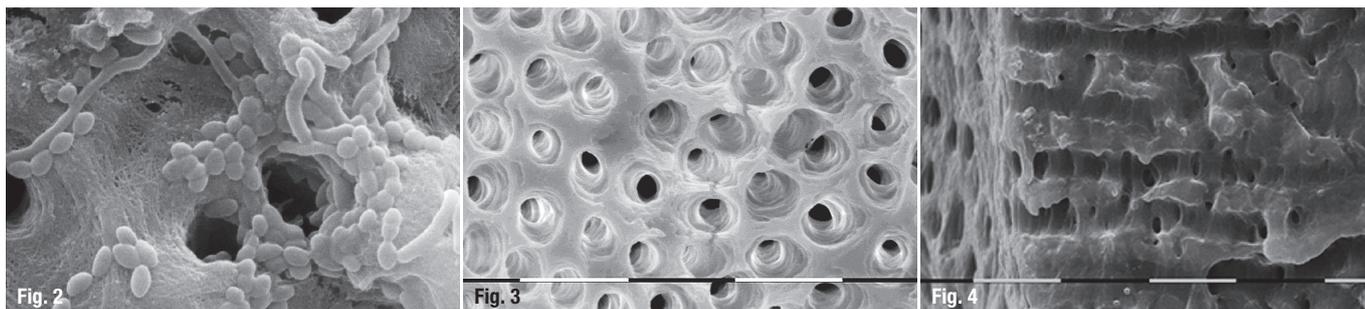

## Laser-enhanced irrigation

Laser energy has been documented to enhance the known effects of NaOCl irrigation through both heating the solution within the canal system and its distant antibacterial effects. But not all laser wavelengths have demonstrated to be equal in effectiveness. The best effects are when NaOCl is combined with an Er:YAG laser as compared to NaOCl alone or when utilised with other laser types.<sup>3</sup> Antibacterial effects were reported to be the best with this combination of irrigant and laser.<sup>4</sup> The higher wavelength of the Er:YAG compared to the Nd:YAG or diode was more effective in smear layer removal, hence, better at bacterial elimination within the canal system.<sup>5</sup> Utilisation of a EDTA as an irrigant alternated with NaOCl provides the best debridement of the canal system with enhancement with a Er:YAG laser, as these two solutions have a synergistic effect complementing each other's effects in the canal anatomy.<sup>6</sup>

Additionally, the Er:YAG laser (LiteTouch™, AMD LASERS) creates hydrodynamic pressure following cavitation bubble expansion and collapse when the irrigation solution is activated in the chamber.<sup>7-9</sup> Placement of the laser tip does not require entry into the canals to achieve the desired effects and activation of the irrigation solution in the chamber is sufficient to affect the entire canal system. The LiteTouch™ Er:YAG laser energy is set at a sub-ablative power level which allows its use without structural changes to the hard tissue within the tooth. This eliminates the risks of ledging and perforation of the pulpal floor allowing safe usage within the tooth.

When the Er:YAG laser is activated, a heat pulse is generated by the laser radiation delivered via a sapphire tip into an absorbing liquid (irrigant). This results in tensile stress with cavitation being induced in the liquid in front of the sapphire tip at a distance far below the optical penetration depth of the laser radiation. Bubble expansion and collapse cause the surrounding fluid to flow at a speed of up to 12m/s travelling throughout the canal system. This causes rapid displacement of intracanal fluid via radial and longitudinal pressures sufficient to drive irrigants into the canal anatomy and clean the dentinal tubules

significantly. This photomechanical activation of the irrigants includes a temperature rise in the irrigant increasing its effectiveness in debridement of dentinal walls and its tubules and increases the chemical properties of the irrigants.

## LiteTouch-Induced Photomechanical Irrigation (LT-IPI)

Endodontic treatment is initiated with access to the pulp chamber, which may be performed by traditional methods using burs or by ablation of the enamel and dentine with the LiteTouch™ Er:YAG laser. As the laser is ineffective in removal of ceramics and metals, such as those used in fixed prosthetics and also amalgam, carbides and diamonds are needed create access through these materials. Once dentine has been reached the laser may be utilised to unroof the pulp chamber (hard tissue mode). An additional benefit of the Er:YAG laser to access the pulp chamber is that it provides decontamination and removal of bacterial debris and pulpal tissue to yield a cleaner chamber aiding it identification of the canal orifices (soft tissue mode).

Once the canal orifices are identified, hand files are utilised to establish a glide path to the apical working length in each canal. Canals are then enlarged to the desired ISO canal size with either hand or rotary files (Fig. 5a). Laser-assisted canal irrigation requires canal preparation to an apical preparation ISO 25/30 at a minimum. A canal taper of .04 or .06 for the final instrumentation is recommended. Sodium hypochlorite (NaOCl) is utilised within the chamber and canals during instrumentation both as a pulpal tissue solvent and to lubricate the files within the canal, de-


**Fig. 2:** SEM showing bacteria and pulpal debris in the apical 1/3 that was not removed fully using standard irrigation protocol.

**Fig. 3:** SEM showing complete removal of bacteria and pulpal tissue in the apical 1/3 after irrigation using the LT-IPI™ protocol.

**Fig. 4:** SEM cross-section showing complete removal of bacteria and pulpal tissue in the apical 1/3 after irrigation using the LT-IPI™ protocol leaving dentine tubules open.

(All pictures © Prof. Georgi Tomov, Plovdiv, Bulgaria)

**Figs. 5a–c:** LiteTouch-Induced Photomechanical Irrigation protocol (LT-IPI).

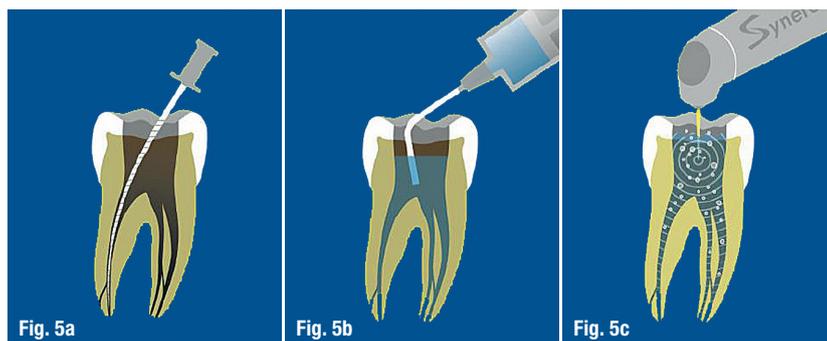


**Fig. 6:** Accessory anatomy evident in the apical that has been filled with sealer accessible due to use of the LiteTouch Er:YAG laser.

(© Dr David Guex, Lyon, France)

**Fig. 7:** Accessory apical anatomy filled with sealer due to use of the LiteTouch™ Er:YAG laser.

(© Prof. Georgi Tomov, Plovdiv, Bulgaria)

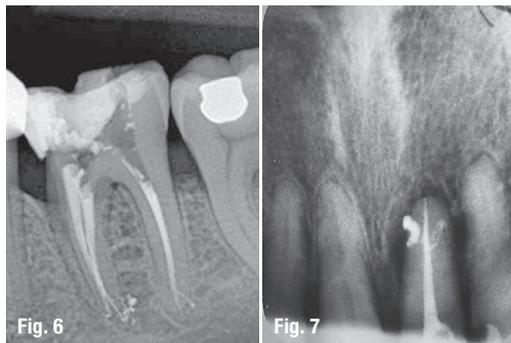

creasing the potential of file separation that can occur when instrumenting a dry canal (Fig. 5b).

Photo-activation of the irrigant within the canal system is performed using the Er:YAG laser with a 0.4/17 or 0.6/17 mm tip which assists in removal of the debris created by the files. Between each rotary file, the chamber is filled with NaOCl and the tip of the laser is placed into the chamber and the solution activated with the laser at 40 mJ at 10 Hz with an average power of only 0.5 W for 20 seconds (Fig. 5c). The chamber is suctioned and fresh NaOCl is placed into the tooth and the next file is used for instrumentation. It is unnecessary to place the lasers tip into the canals themselves, as activation of the solution within the chamber transmits down the irrigant into the canals to the apical aspect of the roots. Laser activation may also be performed with 17% EDTA solution alternated with NaOCl. The benefit of EDTA solution is its chelation effect opening canal anatomy so that the next round of NaOCl can reach more pulpal tissue not accessible to the files in fins, as well as accessory and lateral canals.

Following final instrumentation of the canals with rotary files, the chamber is filled with NaOCl and the Er:YAG tip is placed into the chamber again and activated for a minimum of 60 seconds. This allows the photo-activated irrigant to clear debris and remaining pulpal tissue from the complete canal system. The

irrigation solution is suctioned from the chamber and fresh irrigant placed and photo-activation repeated until no visible debris (cloudiness) is noted in the chamber fluid. This indicated that all accessible debris has been removed from the canal system. Any remaining solution is suctioned from the tooth and the canals are dried with paper points. Obturation is then accomplished using the practitioners preferred method and materials allowing obturation of anatomy inaccessible by instrumentation with files (Figs. 6 & 7).

## Conclusion

The key to endodontic success is two pronged, cleaning the system and sealing it. Although, rotary files have improved the efficiency of instrumentation they are unable to reach any more of the anatomy that handfiles are able to reach. Cleaning of the canal system is keyed to irrigation of the canal system to improve debris removal in anatomy that the files are unable to contact. When anatomy is not fully cleaned, sealer is unable to fill this leaving bacteria to inhabit those areas which may lead to endodontic failure over time. Laser enhanced activation of endodontic irrigants cleans more anatomy adjacent to the main canals so that a more complete obturation of the canal system can occur. An added benefit is that the laser has an antibacterial effect, killing bacteria within the canal anatomy as well as distant to where the irrigation solution may reach essentially sterilising the entire tooth to the periodontal ligament.

## contact

**Dr Gregori M. Kurtzman**  
 Implant Cosmetic Dental Center  
 3801 International Drive 102  
 Silver Spring, MD 20906, USA  
 Tel.: +1 301 5983500  
[dr\\_kurtzman@maryland-implants.com](mailto:dr_kurtzman@maryland-implants.com)  
[www.maryland-implants.com](http://www.maryland-implants.com)

Author details



## Kurz & bündig

Der Schlüssel für eine erfolgreiche Endodontie liegt einerseits in der gründlichen Reinigung des Wurzelkanalsystems und andererseits in dessen Versiegelung. Obwohl rotierende Feilen die Effizienz der Instrumentierung verbessert haben, kann der Behandler feine und verästelte Nebenkanäle mechanisch nur schwer erreichen. Um bakterielle Ablagerungen an schwer zugänglichen Stellen zu beseitigen, erfolgt die Reinigung des Kanalsystems daher mittels Natriumhypochloritlösungen. Doch auch unter Einsatz von Spüllösungen lassen sich nicht alle Bakterien entfernen. Studien haben gezeigt, dass die zusätzliche Anwendung eines Er:YAG-Lasers die Effizienz dieser Lösungen verstärkt und die Desinfektion im Wurzelkanalsystem verbessert. Neben den Hauptkanälen können durch Laseranwendung auch die stark verzweigten Nebenkanäle und andere schwer zugängliche Bereiche gesäubert werden, sodass eine vollständige Obturation des Kanalsystems stattfinden kann. Ein zusätzlicher Vorteil ist die antibakterielle Wirkung des Lasers; hierdurch werden Bakterien innerhalb der Kanalanatomie abgetötet und der gesamte Zahn bis hin zum parodontalen Ligament sterilisiert.

Literature



# SUBSCRIBE NOW!

ceramic implants—international magazine of ceramic implant technology

Interdisciplinary, transparent & market orientated

ONLINE  
SUBSCRIPTION HERE



www.oemus-shop.de

laser 4/17



## stay updated:

research

**From titanium to zirconia implants**

industry

**The DWOS Synergy™ workflow in implantology**

technology

**Shifting of dental implants through ISO standards**

Fax: **+49 341 48474-290**

I would like to subscribe for the following journals:

- ceramic implants—international magazine of ceramic implant technology 2 issues p.a. € 30\*
- implants—international magazine of oral implantology 4 issues p.a. € 44\*

Terms & Conditions: The subscription may be cancelled in written form without due justification within 14 days of order by contacting OEMUS MEDIA AG, Holbeinstraße 29, 04229 Leipzig, Germany. Dispatching notification in good time will suffice. The subscription is automatically extended by another twelve months if it is not cancelled in written form 6 weeks prior to the end of the reference period.

\* All prices include VAT, shipping and handling (within Germany).

\_\_\_\_\_  
Last Name, First Name

\_\_\_\_\_  
Company

\_\_\_\_\_  
Street

\_\_\_\_\_  
ZIP/City/Country

\_\_\_\_\_  
E-Mail

\_\_\_\_\_  
Signature

\_\_\_\_\_  
Credit Card Number

\_\_\_\_\_  
Expiration Date

\_\_\_\_\_  
Security Code